

WIDEX MENU™

THE CHOICE IS YOURS



WIDEX MENU™





IMAGINE IF YOU COULD
OFFER YOUR CLIENTS
A **BASIC HEARING** AID
THAT NOT ONLY SOUNDS
GREAT, BUT CAN ALSO
BE **TAILOR-MADE** TO
SUIT THEIR **HEARING**
LOSS AND **LIFESTYLE**.

MENU, THE NEW
HEARING AID SOLUTION
FROM WIDEX DOES
JUST THAT.



THE **BASIC** CHOICE

The basic MENU™ option is a 3-channel hearing aid that offers more technology than ever before seen at this price.

It includes our unique sound processing platform, noise reduction for cleaner, crisper sound, high level compression for comfortable easy-listening sound and state-of-the art feedback cancellation.

The basic MENU option also provides quick, instant fitting to make it easy for you and your clients.

POPULAR MODELS



Choose between a wide range of outstanding models and stylish colours. MENU is available in BTE (9, 19) and ITE models, as well as the micro *m*. And introducing the new, small and attractive invisible-in-canal model.

ADVANCED WIDEX TECHNOLOGY

The Flex platform – a new, flexible and improved sound platform – features a very fast chip for better performance and better sound.

FAST AND EASY FITTING

Instant fitting with our COMPASS software means you can have MENU up and running in no time. And our instant ear-tips give you extra flexibility.

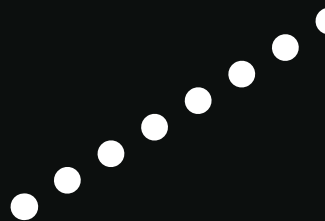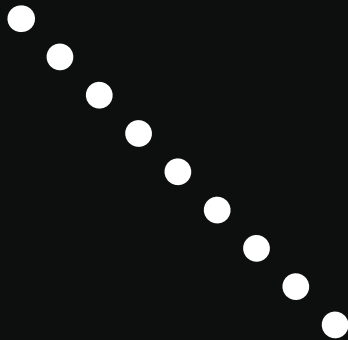

THE **FLEXIBLE** CHOICE

MENU also lets you offer your clients the flexibility of 5 and 10 channels to better suit their needs and preferences.

You can also help them choose from a large range of add-on features such as the HD Locator with Speech Tracer which increases speech understanding, or the Audibility Extender which extends the range of audible sounds.

With WIDEX MENU, the choice is yours.

POPULAR MODELS

Choose between a wide range of outstanding models and stylish colours. MENU is available in BTE (9, 19) and ITE models, as well as the micro *m*. And introducing the new, small and attractive invisible-in-canal model.

TAILOR-MADE

Quickly and easily match the right MENU to your client's lifestyle.

ADVANCED WIDEX TECHNOLOGY

The Flex platform – a new, flexible and improved sound platform – features a very fast chip for better performance and better sound.

USER-FRIENDLY ACCESSORIES

From mobile phone connection to remote control, users can enjoy a wide range of accessories.

FAST AND EASY FITTING

Instant fitting with our COMPASS software means you can have MENU up and running in no time. And our instant ear-tips give you extra flexibility.

FLEXIBILITY

Your clients only pay for what they need.



THE CHOICE IS IN THE FEATURES

With MENU add-ons, you can help users choose unique Widex features to match their hearing needs and their budget.



HD LOCATOR WITH SPEECH TRACER

This directional microphone system aids noise reduction and helps users understand speech from varying directions; the HD Locator is fully adaptive so it produces the best signal-to-noise ratio for the relevant listening environment.



AUDIBILITY EXTENDER

The Audibility Extender expands the range of audible sounds by moving inaudible sounds – like high-frequency speech and environmental sound – down to lower frequencies where they can be heard.



SPEECH ENHANCER

The Speech Enhancer dampens unwanted noise and increases speech intelligibility – all while providing a comfortable listening experience.



ZEN

The popular Zen music and tone program not only helps users to relax but has been proven to help with the effects of tinnitus.



SOUND DIARY

The Sound Diary allows you to easily monitor users' data and hearing aid use – including characteristics of their listening environment and program selection – so you can fine-tune and provide better counselling.



T-DEX

An affordable and easy solution that connects hearing aids with mobile phones and provides excellent sound; connects to MENU via a telecoil.



SMARTSPEAK

SmartSpeak uses actual recorded speech to inform users of various hearing aid functions such as low battery and program selection.



EXTRA PROGRAMS

A range of extra programs can be added to suit the user's listening situation.



THE CHOICES FOR EVERY LIFESTYLE

Your clients all have varying lifestyles, and the different MENU combinations make it easy for you to find the perfect solution for them.

Whether they are at work, active playing sport, relaxing with friends or just enjoying the comforts of home – your clients can enjoy MENU in lots of different situations. You can simply customise MENU for each of your clients as they prefer.

The choice is up to you.







By choosing Widex hearing aids, you are choosing products from a company that has been WindMade certified. WindMade is the first global consumer label identifying companies that use wind power.

Widex, COMPASS and MENU are trademarks of Widex A/S.



Printed by KLS/2012-02
9 502 2537 001 #01

Widex A/S, Nymoellevej 6, DK-3540 Lyngø

WWW.WIDEX.PRO

WIDEX[®]
HIGH DEFINITION HEARING

FEATURE OVERVIEW

Whether you offer your client the basic MENU, or provide them with a range of add-ons – the choice is yours.

	BASIC MENU	ADD-ON
LISTENING IN QUIET		
High Level Compression (WDRC/EDRC available)	●	
Automatic Output Control	●	
Sound Stabilisers	●	
LISTENING IN NOISE		
Noise Reduction SIS	●	
Omni/fixed directional microphones	●	
HD Locator with Speech Tracer		○
Speech Enhancer		○
ENHANCED USER COMFORT		
Multidirectional Active Feedback Cancelling	●	
AISA	●	
EXTENDED AUDIBILITY		
Audibility Extender		○
USER CONTROLS		
SmartTone	●	
Listening programs	2	
Extra programs		○
Zen*		○
SmartSpeak		○
RC4-1		○
RC4-2		○
T-DEX**		○
FITTING AND FINE TUNING		
Basic Sensogram	●	
Expanded Sensogram	●	
In-situ vent effect	●	
Occlusion manager	●	
Sound Diary/Event log		○
PROGRAMMING		
NOAHlink	●	
iP5	●	
Hi-PRO	●	



● YOUR BASIC MENU

○ ADD-ONS: FEATURES THAT CAN BE ADDED

*Zen is not available with the iP5 programmer **T-DEX is available in certain countries only

MENU ADD-ONS

Choosing a feature to suit your client could not be easier.



HD LOCATOR WITH SPEECH TRACER

This directional microphone system aids noise reduction and helps users understand speech from varying directions; the HD Locator is fully adaptive so it produces the best signal-to-noise ratio for the relevant listening environment.



SPEECH ENHANCER

The Speech Enhancer dampens unwanted noise and increases speech intelligibility – all while providing a comfortable listening experience.



SOUND DIARY

The Sound Diary allows you to easily monitor users' data and hearing aid use – including characteristics of their listening environment and program selection – so you can fine-tune and provide better counselling.



SMARTSPEAK

SmartSpeak uses actual recorded speech to inform users of various hearing aid functions such as low battery and program selection.



AUDIBILITY EXTENDER

The Audibility Extender expands the range of audible sounds by moving inaudible sounds – like high-frequency speech and environmental sound – down to lower frequencies where they can be heard.



ZEN

The popular Zen music and tone program not only helps users to relax but has been proven to help with the effects of tinnitus.



T-DEX

An affordable and easy solution for connecting hearing aids with mobile phones; connects to MENU via a telecoil.



EXTRA PROGRAMS

A range of extra programs can be added to suit the user's listening situation.

COLOURS AND MODELS:



Warm beige (021)



Jet black (030)



Cappuccino brown (071)



Winter silver (072)



Titan grey (073)



Charcoal grey (035)



Light beige (025)



Dark brown (027)









Clay brown (028)



Tan silk (081)

FITTING RANGE:

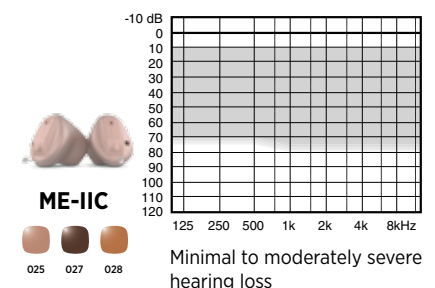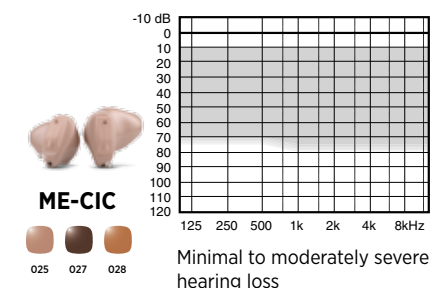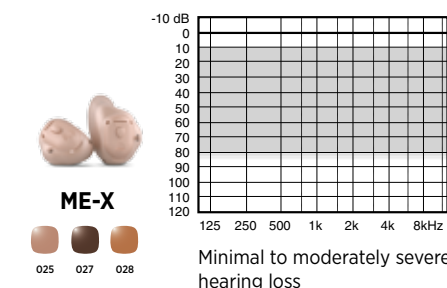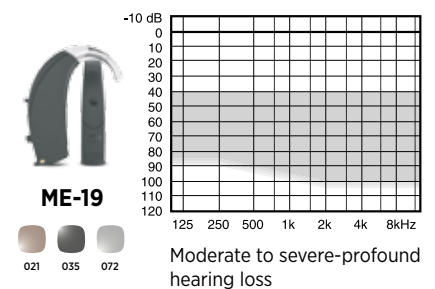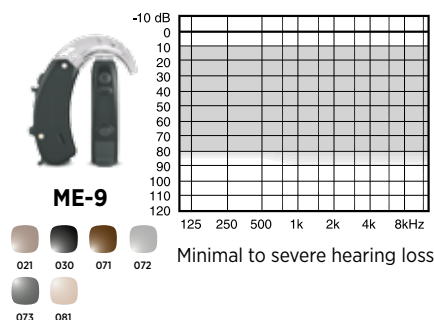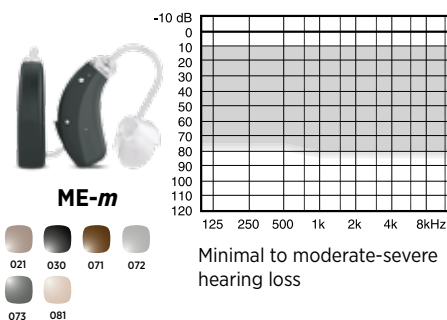

INSTANT EAR-TIP OPTIONS FOR *m* AND 9 MODELS



Instant open ear-tip



Instant tulip ear-tip



Instant double ear-tip

Not shown in real size